


Adaptor for slides for macro-sections

General features

PEPT (polyethylene terephthalate) adaptor for the staining of big slides for macro-sections with standard slide boxes.

Code	Packaging
20-E103/SL/R30A	2 pcs.

Manufacturer: Bio-Optica Milano S.p.A.

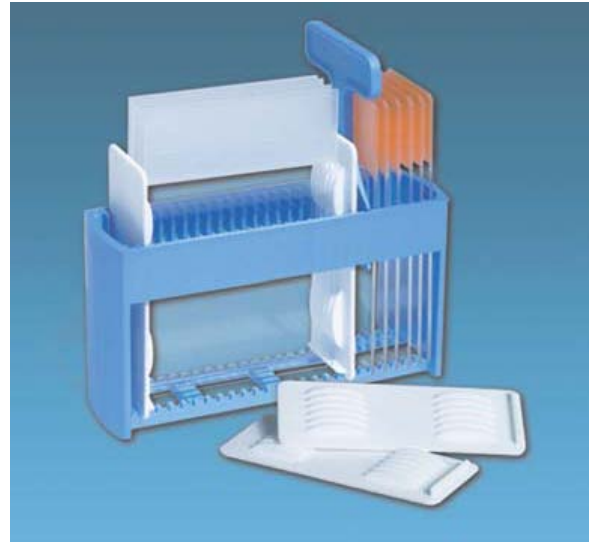

Technical features

Plastic slides' adaptor, resistant to staining solutions, solvents and reagents generally used in histology labs (like alcohols and xylene).

The adaptors have to be used in couple (like in the picture) and each couple allows the insert of 5 slides for macro-sections.

Each adaptor has this dimensions: 27 x 67 mm and thickness 1,3 mm to be well-adjusted in the standard slide boxes. In this way, it is possible to stain together in the same slide box standard slides and big slides.

Packaging: 2 pcs. in a cellophane box enough to guarantee the quality during the time.

For this kind of product, particular rules for the transport and the storage are not necessities.

Expiry date: not applicable.